

### **Related Topics**

Maxwell's equations, wire loop, flat coils, Biot-Savart's law, Hall effect.

#### Principle

The spatial distribution of the field strength between a pair of coils in the Helmholtz arrangement is measured. The spacing at which a uniform magnetic field is produced is investigated and the superposition of the two individual fields to form the combined field of the pair of coils is demonstrated.

#### Equipment

1 Pair of Helmholtz coils	06960.00
1 Power supply, universal	13500.93
1 Digital multimeter	07134.00
1 Teslameter, digital	13610.93
1 Hall probe, axial	13610.01
2 Meter scale, demo, I = 1000 mm	03001.00
1 Barrel base -PASS-	02006.55
1 Support rod -PASS-, square, I = 250 mm	02025.55
1 Right angle clamp -PASS-	02040.55
3 G-clamp	02014.00
1 Connecting cord, I = 750 mm, blue	07362.04
3 Connecting cord, I = 750 mm, red	07362.01



Fig. 1: Set-up of experiment P2430301

TESS expert PHYWE

#### Tasks

- 1. Measure the magnetic flux density along the *z*-axis of the flat coils when the distance between them a = R (R = radius of the coils) and when it is greater and less than this.
- 2. Measure the spatial distribution of the magnetic flux density when the distance between coils a = R, using the rotational symmetry of the set-up:
  - a. measurement of the axial component  $B_z$
  - b. measurement of radial component B<sub>r</sub>
- 3. Measure the radial components  $B_r^{*}$  and  $B_r^{"}$  of the two individual coils in the plane midway between them and to demonstrate the overlapping of the two fields at  $B_r = 0$ .

#### Set-up and Procedure

Connect the coils in series and in the same direction, see Fig. 2; the current must not exceed 3.5 A (operate the power supply as a constant current source). Measure the flux density with the axial Hall probe (measures the component in the direction of the probe stem).

The magnetic field of the coil arrangement is rotationally symmetrical about the axis of the coils, which is chosen as the *z*-axis of a system of cylindrical coordinates (z, r,  $\Phi$ ). The origin is at the centre of the system. The magnetic flux density does not depend on the angle  $\Phi$ , so only the components  $B_z$  (z, r) and  $B_r$ (z, r) are measured.

Clamp the Hall probe on to a support rod with barrel base, level with the axis of the coils. Secure two rules to the bench (parallel or perpendicular to one another,


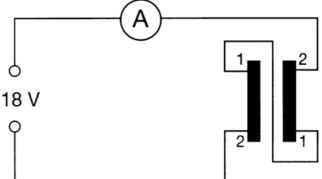


Fig. 2: Wiring diagram for Helmholtz coils.



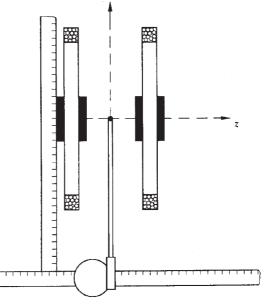

see Figs. 3–5). The spatial distribution of the magnetic field can be measured by pushing the barrel base along one of the rules or the coils along the other one.


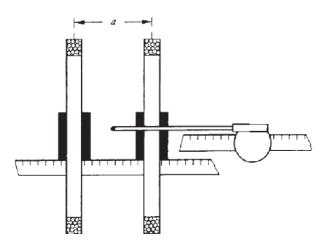


Fig. 3: Measuring B(z, r = 0) at different distances *a* between the coils.

**PHYWE** excellence in science



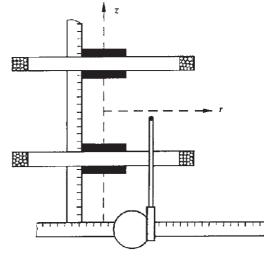

## Notes

Always push the barrel base bearing the Hall probe along the rule in the same direction.

- 1. Along the *z*-axis, for reasons of symmetry, the magnetic flux density has only the axial component  $B_z$ . Fig. 3 shows how to set up the coils, probe and rules. (The edge of the bench can be used instead of the lower rule if required.) Measure the relationship B(z, r = 0) when the distance between the coils a = R and, for example, for a = R/2 and a = 2R.
- 2. When distance a = R the coils can be joined together with the spacers. a) Measure  $B_z(z, r)$  as shown in Fig. 4. Set the *r*-coordinate by moving the probe and the *z*-coordinate by moving the coils. Check: the flux density must have its maximum value at point (z = 0, r = 0). b) Turn the pair of coils through 90° (Fig. 5). Check the probe: in the plane  $z = 0, B_z$  must = 0.
- 3. Short-circuit first one coil, then the other. Measure the radial components of the individual fields at z = 0.

# Theory and evaluation

From Maxwell's equation

$$\oint_{K} \vec{H} d\vec{s} = I + \int_{F} \int \vec{D} d\vec{f} dt$$
(1)

where K is a closed curve around area F, we obtain for direct currents ( $\dot{D}$  = 0), the magnetic flux law

$$\oint_{K} \vec{H} d\vec{s} = I$$
<sup>(2)</sup>

which is often written for practical purposes in the form of Biot-Savart's law:

$$d\vec{H} = \frac{I}{4\pi} \frac{d\vec{\iota} \times \vec{\rho}}{\rho^3} \tag{3}$$

where  $\vec{\rho}$  is the vector from the conductor element  $d\vec{i}$  to the measurement point and  $d\vec{H}$  is perpendicular to both these vectors.

The field strength along the axis of a circular conductor can be calculated using equation (3). (Fig. 6).

The vector  $d\vec{i}$  is perpendicular to, and  $\vec{\rho}$  and  $d\vec{H}$  lie in, theplane of the sketch, so that

$$dH = \frac{I}{4\pi\rho^3} d\iota = \frac{I}{4\pi} \cdot \frac{d\iota}{R^2 + z^2}$$
(4)

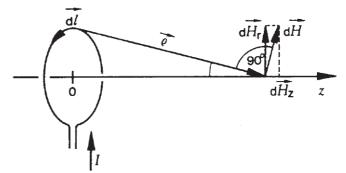


Fig. 6: Sketch to aid calculation of the field strength along the axis of a wire loop.

TEP .3.03	Magnetic field of paired coils	TESS PHYWE
-01	in Helmholtz arrangement	expert

 $d\vec{H}$  can be resolved into a radial  $dH_r$  and an axial  $dH_z$  component. The  $dH_z$  components have the same direction for all conductor elements and the quantities are added; the  $dH_r$  components cancel one another out, in pairs. Therefore,

$$H_r = 0 \tag{5}$$

and

$$H = H_z = \frac{I}{2} \cdot \frac{R^2}{(R^2 + z^2)^{3/2}}$$
(6)

along the axis of the wire loop, while the magnetic flux density

$$B(z) = \frac{\mu_0 \cdot I}{2R} \cdot \frac{1}{\left(1 + \left(\frac{z}{R}\right)^2\right)^{3/2}}$$
(7)

The magnetic field of a flat coil is obtained by multiplying (6) by the number of turns *N*. Therefore, the magnetic flux density along the axis of two identical coils at a distance  $\alpha$  apart is

$$B(z,r=0) = \frac{\mu_0 \cdot I \cdot N}{2R} \cdot \left(\frac{1}{\left(1 + {A_1}^2\right)^{3/2}} + \frac{1}{\left(1 + {A_2}^2\right)^{3/2}}\right)$$
(8)

where

$$A_1 = \frac{z + \alpha/2}{R}, \ A_2 = \frac{z - \alpha/2}{R}$$

When z = 0, flux density has a maximum value when  $\alpha < R$  and a minimum value when  $\alpha > R$ . The curves plotted from our measurements also show this (Fig. 7); when  $\alpha = R$ , the field is virtually uniform in the range

$$-\frac{R}{2} < z < +\frac{R}{2}$$

Magnetic flux density at the mid-point when  $\alpha = R$ :

$$B(0.0) = \frac{\mu_0 \cdot I}{2R} \cdot N \cdot \frac{2}{\left(\frac{5}{4}\right)^{\frac{3}{2}}} = 0.716 \,\mu_0 \cdot N \cdot \frac{I}{R}$$

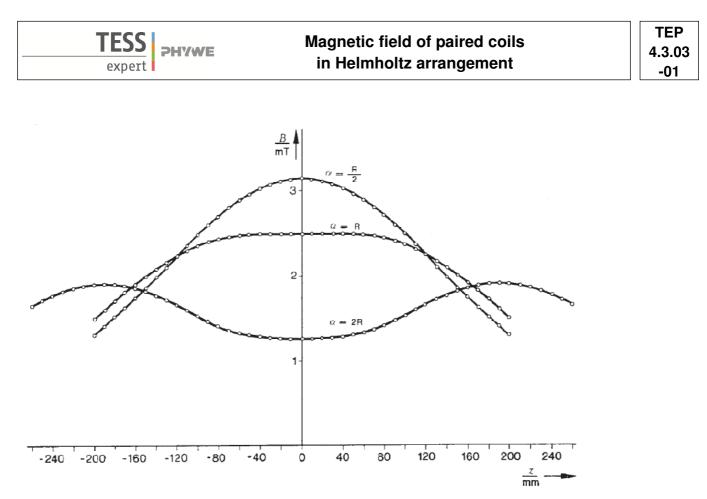


Fig. 7: B(r = 0) as a function of z with the parameter  $\alpha$ .

when N = 154, R = 0.20 m and I = 3.5 A this gives:

B(0.0) = 2.42 mT.

Our measurements gave B(0.0) = 2.49 mT.

Figs. 8 and 9 shows the curves  $B_z(z)$  and  $B_r(z)$  measured using *r* as the parameter; Fig. 10 shows the super-position of the fields of the two coils at  $B_r = 0$  in the centre plane z = 0.

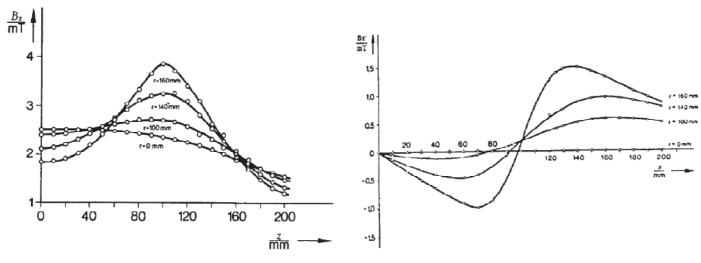
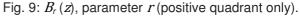


Fig. 8:  $B_z(z)$ , parameter r (positive quadrant only).



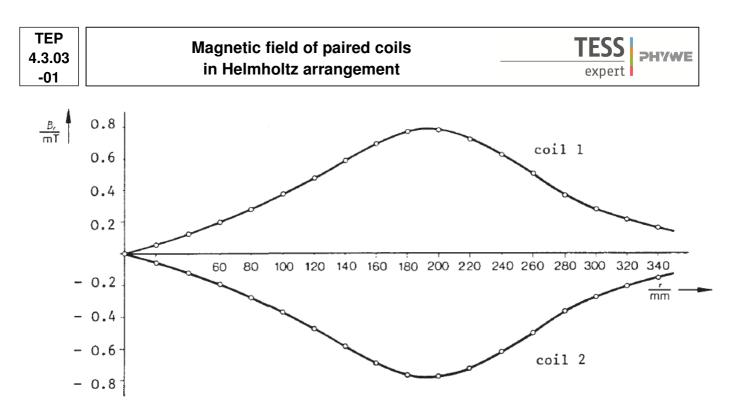


Fig. 10: Radial components  $B_{r'}(r)$  and  $B_{r''}(r)$  of the two coils when z = 0.