

Malus' law (Item No.: P2250400)

Curricular Relevance



Difficulty



Difficult

Preparation Time



10 Minutes

Execution Time



20 Minutes

Recommended Group Size



3 Students

Additional Requirements:

Experiment Variations:

Keywords:

Electric theory of light, polarization, polarizer, analyzer, Brewster's law, Malus' law

Introduction

Overview

Linear polarized light passes through a polarization filter. Transmitted light intensity is determined as a function of the angular position of the polarization filter.



Fig. 1: Experimental set-up: Malus' law.

Student's Sheet

Printed: 10.08.2017 12:01:46 | P2250400

Equipment

Position No.	Material	Order No.	Quantity
1	Laser, He-Ne 1.0 mW, 220 V AC	08181-93	1
2	Optical profile bench, l = 60 cm	08283-00	1
3	Base for optical profile bench, adjustable	08284-00	2
4	Slide mount f. opt. pr.-bench, h = 30 mm	08286-01	3
5	Polarising filter on stem	08610-00	1
6	Photoelement f. opt. base plate	08734-00	1
7	Digital multimeter	07122-00	1

Tasks

1. The plane of polarization of a linear polarized laser beam is to be determined.
2. The intensity of the light transmitted by the polarization filter is to be determined as a function of the angular position of the filter.
3. Malus' law must be verified.

Set-up and procedure

The experiment is set up according to Fig. 1. It must be made sure that the photocell is totally illuminated when the polarization filter is set up.

If the experiment is carried out in a non darkened room, the disturbing background current i_0 must be determined with the laser switched off and this must be taken into account during evaluation.

The laser should be allowed to warm up for about 30 minutes to prevent disturbing intensity fluctuations.

The polarization filter is then rotated in steps of 5° between the filter positions $\pm 90^\circ$ and the corresponding photo cell current (most sensitive direct current range of the digital multimeter) is determined.

Theory and evaluation

Let AA' be the Polarization planes of the analyzer in Fig. 2. If linearly polarized light, the vibrating plane of which forms an angle φ with the polarization plane of the filter, impinges on the analyzer, only the part

$$E_A = E_0 \cdot \cos\varphi \tag{1}$$

will be transmitted.

As the intensity I of the light wave is proportional to the square of electric field intensity vector \vec{E} , the following relation (Malus' law) is obtained:

$$I_A = I_0 \cdot \cos^2 \varphi \tag{2}$$


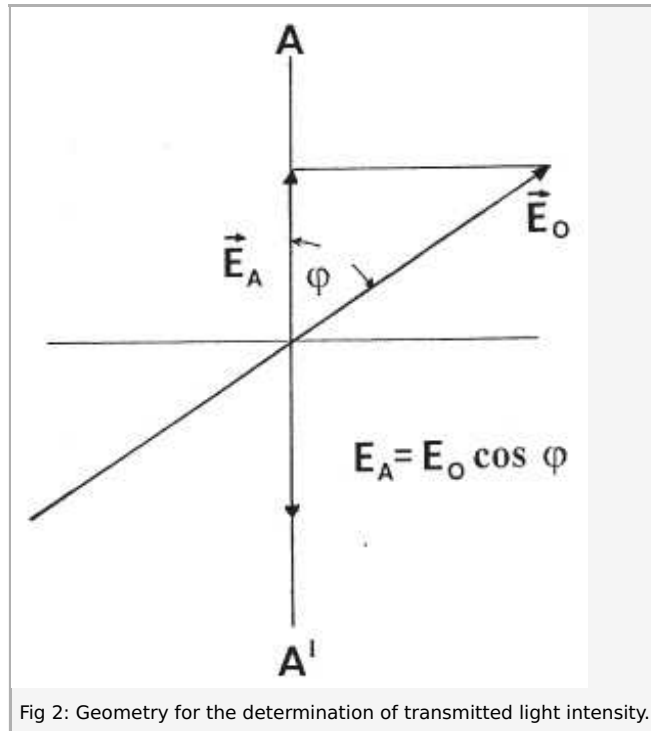


Fig. 3 shows the photo cell current after background correction (this is a measure for the transmitted light intensity) as a function of the angular position of the polarization plane of the analyzer. The intensity peak for $\varphi = 50^\circ$ shows that the polarization plane of the emitted laser beam has already been rotated by this angle against the vertical.

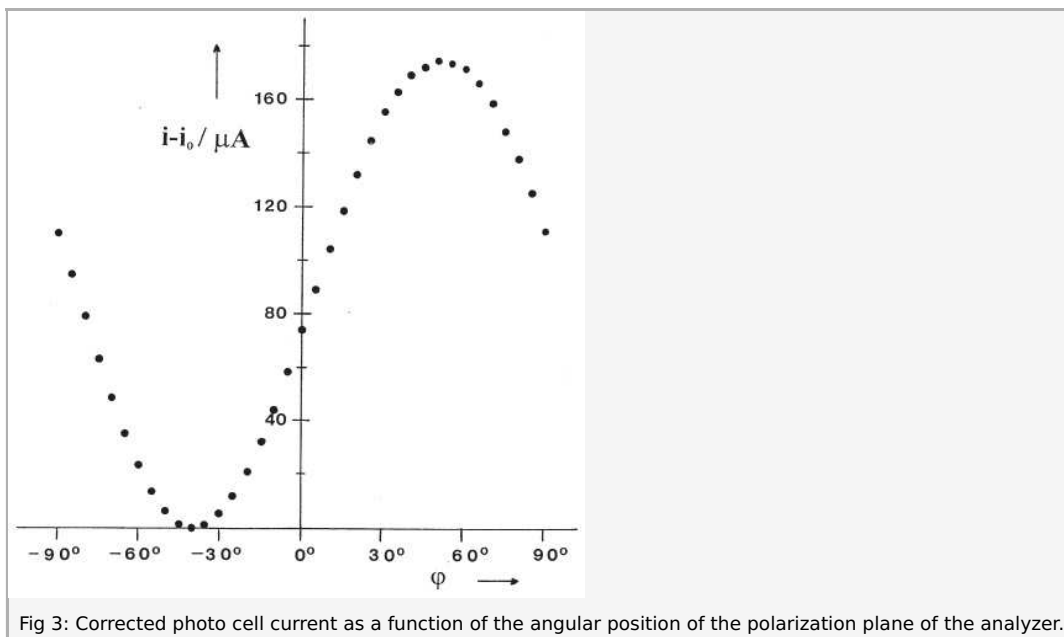

Student's Sheet

Printed: 10.08.2017 12:01:46 | P2250400

Fig. 4 shows the normalized and corrected photo cell currents as a function of the angular position of the analyzer. Malus's law is verified by the initial line's 45° slope (Note: to determine Malus' line in Fig. 4, an angular setting of 50° of the analyzer must be considered for $\varphi = 0^\circ$).

